


DustMatic-2000/2001, Dust Drawer

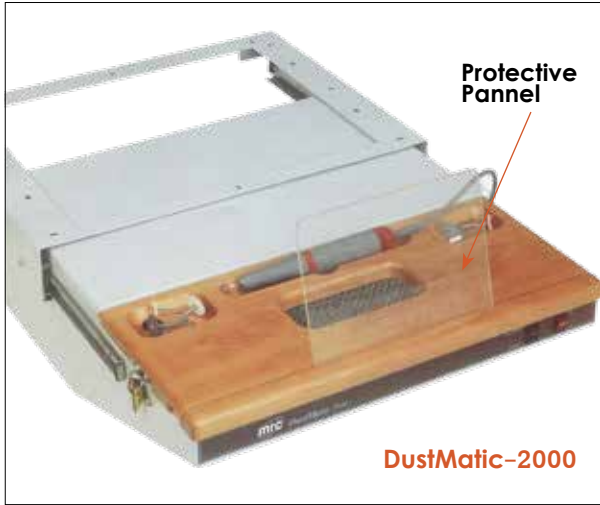

A dust-suction drawer of optimal quality and design, featuring innumerable advantages & improvements in almost all fields of operation, efficiency, user-friendliness, etc.

- Operation is fully automatic: on starting the micromotor, the drawer starts operating, on stopping the micromotor, suction switches off after 20 seconds have elapsed.
- Ergonomic design and construction with first quality materials.
- The working surface is made of beech-wood, giving it a modern look, and allowing for comfortable working.
- Silent, powerful, maintenance-free low pressure induction motor.
- Large-dimensioned, easy-to-change, filter bag (according to international Standards).
- Small-dimensioned drawer, capable of fitting most lab workbenches, with ample leg-room.
- Low operational costs. Easy to install, requires no special tools or training.

Applications:

- Dental
- Ear Phones manufacturer
- Artificial Eye manufacturer
- Archaeology Laboratory
- Beauty manicure
- Jewelry manufacturer.

Include:

Suction drawer complete, with micro-respirable dust filter bag & protective panel.

- Effective suction with large suction bag.
- Silent, maintenance free & applicable for continuous operation.
- Easy mountable under every table.



Model	Power	Air Volume	Size (mm)	Weight	Operation
2000	100W	350m³h	490x510x180	14kg	Auto/Manual
2001	100W	350m³h	490x510x180	13.5kg	Manual only

